

RA640 BASIC

FRONT PROTECTION APRON

NEW

Available as ...

Premium Line - The Full Range
Comfort Line - The Economical Alternative

The basic RA640 model has proved effective for years and is used as a front protection usually for short and rare wearing times in conventional X-ray situations without fluoroscopy.

- High-quality, flexible standard lead material with lead equivalent values of 0.25 mm Pb, 0.35 mm Pb or 0.50 mm Pb
- Ideal fit, maximum comfort for the wearer
- Ideal for short and infrequent wear
- Shoulder pads and preparation for adaptable embroidery as standard equipment
- Type testing and certification in accordance with (EU) 2016/425, IEC 61331-1:2014 / DIN EN 61331-1:2016 and IEC 61331-3:2014 / DIN EN 61331-3:2016



Sizes*, recommended for:

| | | |
|-----|----------------------------------|--------------------------------------|
| XS | Chest girth ≤ 84 cm | Waist/hip girth ≤ 91 cm |
| S | Chest girth approx. 85 - 95 cm | Waist/hip girth approx. 92 - 102 cm |
| M | Chest girth approx. 96 - 106 cm | Waist/hip girth approx. 103 - 112 cm |
| L | Chest girth approx. 107 - 117 cm | Waist/hip girth approx. 113 - 123 cm |
| XL | Chest girth approx. 118 - 128 cm | Waist/hip girth approx. 124 - 134 cm |
| XXL | Chest girth approx. 129 - 139 cm | Waist/hip girth approx. 135 - 145 cm |
| 3XL | Chest girth approx. 140 - 150 cm | Waist/hip girth approx. 146 - 156 cm |


* The size depends on the girth of the chest, waist, or hip, depending on which of the three measurements is the largest.

FRONT PROTECTION APRON


BASIC RA640

Different closure types – You can choose ...


We offer different closure systems for this series of aprons. Choose your favourite:

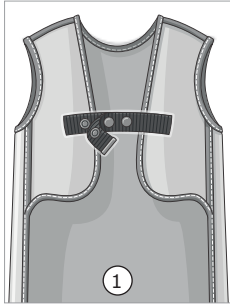


Snap fastener closure
Elastic band with snap fasteners for closing at the back.

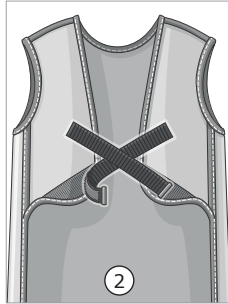

Crossband closure
Crossbands with hook-and-loop fastening at the back.

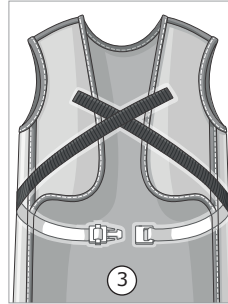

Strap closure
Straps with snap buckle closure for closing at the front.

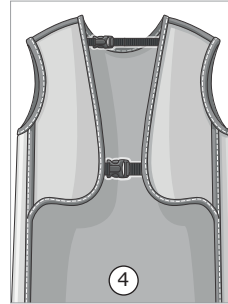

Snap buckle closure
Two webbing straps with snap buckle closures for closing at the back.



Back wings

For model RA640 it is also possible to choose whether protection of 0.25 mm Pb should be maintained in the back wings or dispensed with to reduce the weight of the apron.



Protection/lead equivalent

Novalite, Leadfree or Standard lead



Front, optional:
0.25 mm Pb, 0.35 mm Pb or 0.50 mm Pb

Back, optional:
0.25 mm Pb or without protection

Outer material

Colours (see front page)



Optional:
ComforTex® HPMF, ComforTex® HPMF Hybrid


Sizes

from XS to 3XL (see front page)



Skirt length / Body height:
90 cm / 150 - 160 cm
100 cm / 160 - 170 cm
110 cm / 170 - 180 cm
120 cm / 175 - 190 cm
130 cm / 185 - 195 cm
140 cm / 195 - 205 cm

Optional

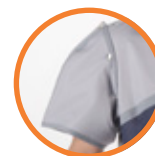

On request:
Sleeves
Breast pocket
Embroidery
Made-to-measure protection

Accessories/Replacement

Bolero
Thyroid/sternum protection



Crossband closure
(Adjustment via snap fasteners)