





RP689

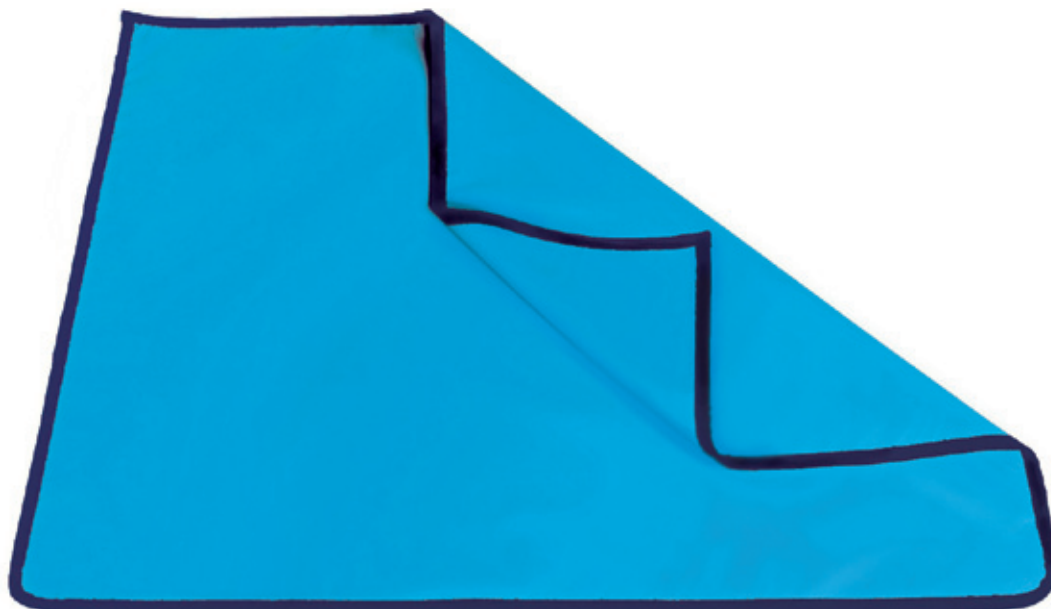

TIP

Do not oversize the dimensions of your patient protective cover. Remember that you have to move and carry the protective cover and that the size stands in direct correlation to the weight.

Radiation protection that is effective for your patient should shield those body parts that are directly adjacent to the area being diagnosed. You may be better equipped with several smaller covers that can be combined depending on the situation.

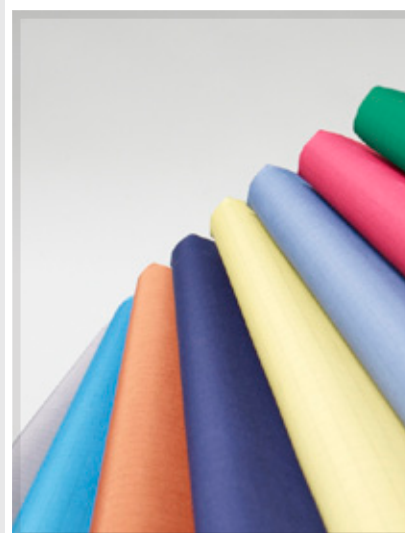

Article No.: **RP689CC*** Patient Cover RP689 (Multi-layer lead rubber material)

Made in Germany



RP689 – Patient Cover

This versatile patient protective cover is a must in all radiology practices: specific protection for your patients, combined with simple handling.

The size can be ordered, depending on your purpose or choice. There are practically no limits. It is also possible to supply the protective cover with hook & burr fasteners or customize to suit your individual needs.

Our MAVIG patient protective cover RP689 is made from high quality, flexible multi-layer lead rubber material with the best radiation protection properties.

- **Protective Material** Multi-layer lead rubber protective material
- **Lead Equivalent** Protection in Pb 0.50 mm or Pb 1.00 mm
- **Outer Cover/Colour** ComforTex® HPMF, available in CC*
- **Size** Our standard widths are 60 cm and 90 cm in your defined length. Specially tailored sizes are also available upon request.
- **Edging** All edges covered with seam band

MAVIG GmbH
X-Ray Protection and
Medical Suspension Systems

Stahlgruberring 5
81829 Munich
Germany

PO Box 82 03 62
81803 Munich
Germany

Phone +49 (0) 89 / 420 96 - 0
Fax +49 (0) 89 / 420 96 - 200
e-Mail info@mavig.com

(*CC - ColourCode: 1 = Curacao, 2 = Orchidee, 3 = Regatta, 4 = Indian Summer, 5 = Limette, 6 = Teddy, 7 = Ocean, 8 = Lagoon)