



EN LED6 MC

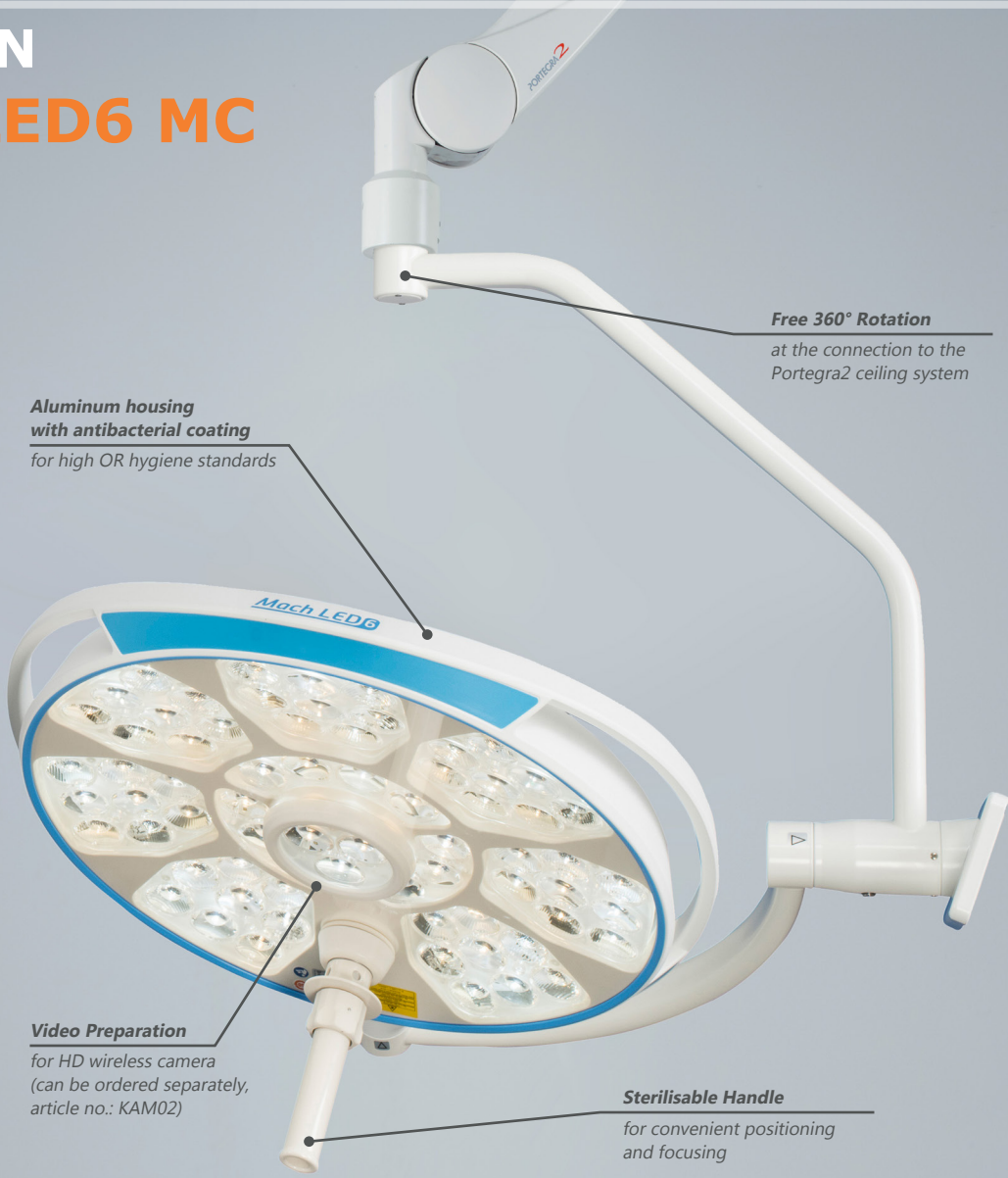

HD wireless camera (optional)

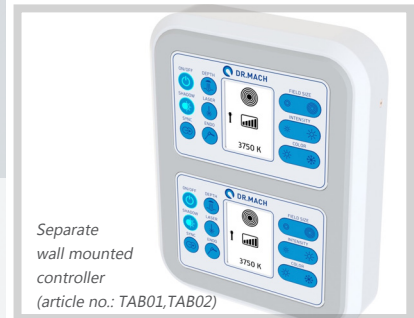

Separate wall mounted controller (article no.: TAB01, TAB02)

- Superior colour rendition
- Enhanced illumination depth
- Large focus field
- Antibacterial coating
- Optimal laminar flow properties
- Laser pointer for correct positioning (optional)
- Camera ready (optional)
- Automatic shadow management (optional)

Models:	LE7069111	LED6 MC, Multi-Colour LED lamp, sterilisable handle and electric preset Portegra2 suspension arm (LE70U628): extension arm 750 mm, external power supply unit
	LE7069111L	LED6 MC, comparable with LE7069111, incl. a laser pointer upgrade
	LE7069111KVS	LED6 MC, comparable with LE7069111, incl. camera preparation and shadow management

Made in Germany **CE, NRTL**

LED6 MC – LED OR Lamp

The LED6 MC is an operating theater light for diagnosis and surgery. The electronic brightness is controlled on the lamp head and the size of the light field can be easily adjusted using the sterilisable handle. Using the optionally available wall panel, the lamp can also be operated remotely. The LED6 MC is available in 3 versions: 1. Basic version, 2. including laser pointer (L) or 3. including camera preparation and shadow management (KVS).

- | | |
|---|--|
| <ul style="list-style-type: none"> • Central illuminance (at 1 m distance) • Colour rendering index R_a at 4500 Kelvin (CRI) • Focusable illuminated field size (at 1 m distance) • Electronic brightness control • Temperature increase in the head area • Low power consumption (total) / Number of LEDs | <p>160,000 lx
 98
 180 - 300 mm / 7.1 - 11.8 in
 30 - 100 % (incl. 5% endoscopy light function)
 0.5 K
 max. 92 VA / 69 LED (66 LED with camera)
 (mains voltage 100 - 240 VAC, 50 - 60 Hz)
 ≥ 60,000 h
 60 - 150 cm / 23.6 - 59.0 in
 117 cm / 46.0 in
 58 cm / 22.8 in
 max. 16.5 kg / 36.4 lbs, depending on model
 IP54 (LE7069111, LE7069111L) / IP53 (LE7069111KVS)
 in accordance with MDR Regulation / CE, NRTL</p> |
| <ul style="list-style-type: none"> • Long service life of lamp • Working range • Height adjustment (at Portegra2 spring arm) • Diameter of lamp head • Weight LED6 MC • Protection class • Medical Product - (EU) 2017/745,
Usage in accordance with DIN VDE 0100-710 / Approval | |

*MC Multi-Colour: Colour temperatures (Kelvin) 3750, 4000, 4250, 4500, 4750

MAVIG GmbH
X-Ray Protection and
Medical Suspension Systems

Stahlgruberring 5
81829 Munich
Germany

PO Box 82 03 62
81803 Munich
Germany

Phone +49 (0) 89 / 420 96 - 0
Fax +49 (0) 89 / 420 96 - 200
e-Mail info@mavig.com