



Victoria RA636

ALL-ROUND PROTECTION: COAT ESPECIALLY FOR WOMEN

The simple, convenient way of putting it on adds to the comfort. The coat is simply closed with a hook-and-loop fastener and sits correctly right away.

Protection/Lead Equivalent Standard Collection

NovaLite, Leadfree
or Standard lead



Front:
0.25 mm Pb, 0.35 mm Pb
or 0.50 mm Pb

Back generally:
0.25 mm Pb

Protection/Lead Equivalent Levio Edition

NovaLite, Leadfree
or Standard lead



Total value Levio-Panel:
0.42 mm Pb or 0.50 mm Pb
(outside 0.25 mm Pb)

Back generally:
0.25 mm Pb

- ✓ X-ray protective coat especially for women
- ✓ Ideal fit, maximum comfort for the wearer
- ✓ Does not require a full overlap at the front
- ✓ Continuous support thanks to the elastic stretch belt (width 10 cm, included)
- ✓ Standard features include shoulder pads
- ✓ Type testing and certification in accordance with (EU) 2016/425, IEC 61331-1:2014 / DIN EN 61331-1:2016 and IEC 61331-3:2014 / DIN EN 61331-3:2016

Alternatively, the Victoria coat can be combined with our shoulder relief system **BRYGGA**.



Sizes*, recommended for:

XS	Chest girth ≤ 84 cm	Waist/hip girth ≤ 98 cm
S	Chest girth approx. 85 - 95 cm	Waist/hip girth approx. 99 - 108 cm
M	Chest girth approx. 96 - 106 cm	Waist/hip girth approx. 109 - 119 cm
L	Chest girth approx. 107 - 117 cm	Waist/hip girth approx. 120 - 130 cm
XL	Chest girth approx. 118 - 128 cm	Waist/hip girth approx. 131 - 141 cm
XXL	Chest girth approx. 129 - 139 cm	Waist/hip girth approx. 142 - 152 cm
3XL	Chest girth approx. 140 - 150 cm	Waist/hip girth approx. 153 - 163 cm

* The size depends on the girth of the chest, waist, or hip, depending on which of the three measurements is the largest.



ALL-ROUND PROTECTION: COAT ESPECIALLY FOR THE FEMALE BODY

VICTORIA RA636

Worldwide Unique!


The Victoria Collection ...

... represents more than just protection - it embodies a modern, sophisticated solution tailored to the needs of the female anatomy. Get ready for the next generation of innovative X-ray protection.

With specially adapted armholes, the Victoria RA636 coat offers ideal protection for the sensitive breast tissue without compromising the freedom of movement.

Princess seams (longitudinal dividing seams for a close-fitting cut) on the front ensure an ideal fit, while a stretch insert in the back ensures an optimal adaptation to the body shape.

Discover the new fusion of functionality and style that combines safety and aesthetic design to meet the needs and demands of women in interventional radiology.

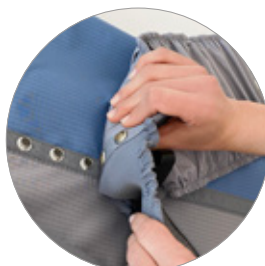

Snap Button Fastening

Optimal adjustment to your waist height is ensured by a variable snap button fastener of the belt.

This allows it to be attached variably in height and to be buttoned to the coat on one or both sides.

Hint:

Please note that sleeves cannot be attached to the Victoria coat due to the adapted armholes. Please combine the Victoria coat with the MAVIG bolero RA637.

Levio®

Optimised Protection Concept



BRYGGA Shoulder Relief



Outlast Temperature Regulation

Higher Protection – Levio Edition

The area with the reinforced protective material is clearly recognisable by its colour.

Optimised Armholes

Ideal protection for sensitive breast tissue

Dividing Seams in the Front

For a good fit but with full radiation protection

Free Movement

For greater freedom of movement, the coat has a cut-out in the back. Of course, the knees are still optimally protected at the front and sides.



Outer Material

Colours
(see front page)



Optional:
ComforTex® HPMF,
ComforTex® HPMF Hybrid

Sizes

from XS to 3XL
(see front page)



Coat length / Body height:
90 cm / 150 - 160 cm
100 cm / 160 - 170 cm
110 cm / 170 - 180 cm
120 cm / 175 - 190 cm
130 cm / 185 - 195 cm
140 cm / 195 - 205 cm

Optional



On request:
Outlast®
BRYGGA®
Breast pocket
Embroidery
Made-to-measure protection

Accessories/Replacement

Bolero
Thyroid/sternum protection



Stretch belt:
(Adjustment via
snap fastener)