



EN

RD639

MAVIG

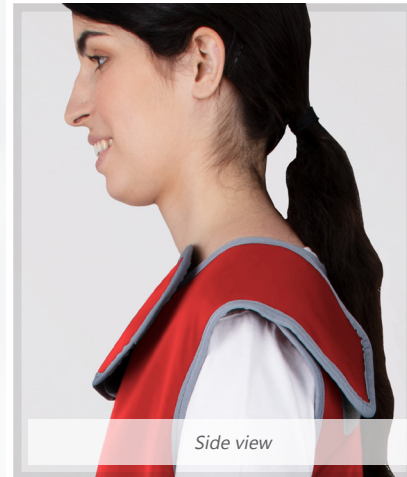

Side view

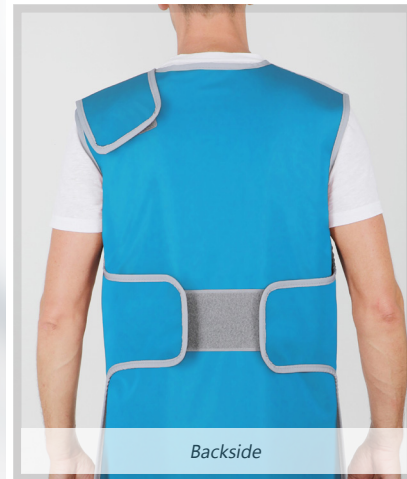

Backside

Model: RD639

Adult dental all-round protection with integrated thyroid protection
Patient protection for digital volume tomography (DVT)

Made in Germany

CE 0302

RD639 – Dental All-Round Protection for digital volume tomography

The RD639 allround protective apron offers optimum radiation protection in digital volume tomography (DVT). The patient is protected at the front as well as the back with the full lead equivalent of 0.35 mm Pb.

The convenient hook and burr fastener on the shoulder allows the apron to be put on comfortably from the side. The circumference can be easily adjusted by means of the hook and burr fasteners in the back. The integrated thyroid protection is height adjustable, so that it does not interfere with the imaging.

Delivery includes special hanger AW106 for one garment.

- **Protective Material** Novalite Plus or Leadfree 2.0
- **Lead Equivalent** 0.35 mm Pb (front and back)
- **Sizes**

Adults:	Small	80.0 cm x 60.0/60.0 cm (L x W, front/back)
	Medium	80.0 cm x 70.0/60.0 cm (L x W, front/back)

(for children please use our children's coat RP668 combined with our thyroid shield RD614)
- **Outer Cover / Colour** ComforTex® HPMF Duplex, available in the colours:
Titan, Blue Admiral, Regatta, Curacao, Arctic, Olive, Lilac, Orchidee, Terra, Onyx, or design: Happy (ComforTex® HPMF only)

MAVIG GmbH

X-Ray Protection and
Medical Suspension Systems

Stahlgruberring 5
81829 Munich
Germany

PO Box 82 03 62
81803 Munich
Germany

Phone +49 (0) 89 / 420 96 - 0
Fax +49 (0) 89 / 420 96 - 200
e-Mail info@mavig.com